



Galaxy Treatment Tables

Introduction

Chattanooga is the world's largest manufacturer of rehabilitation equipment for treating musculoskeletal, neurological and soft tissue disorders. For over six decades, Chattanooga has set a benchmark for leadership, reliability and excellence. Our products contribute to better treatment outcomes in hospitals, clinics and home settings worldwide. We lead the physiotherapy industry by example and, through continuous innovation, provide real world solutions for clinicians and their patients.

As the world's leading brand in rehabilitation, Chattanooga is sold in over 80 countries worldwide. Chattanooga is directly represented in our key markets but also has long lasting relationships with dedicated distributors in the USA, Canada, Latin America, EMEA and Asia Pacific. Our Chattanooga distributors and customers are supported for order placement, training, product delivery and after sales service by a dedicated team of highly skilled DJO employees based all over the world.

Chattanooga is a brand of DJO Global Inc. DJO Global provides solutions for musculoskeletal health, vascular health and pain management. Our products help prevent injuries or rehabilitate after surgery, injury or from degenerative disease, enabling patients to regain or maintain their natural motion.

Galaxy Tables are manufactured by Meden-Inmed Sp. z o.o. exclusively for Chattanooga. A DJO Global Company.

Visit www.DJOGlobal.eu



Contents

Introduction	1
Galaxy Tables	3
Galaxy 5 Section High-Low	4
Galaxy 4 Section High-Low	5
Galaxy 3 Section High-Low	6
Galaxy 2 Section High-Low Head Section	7
Galaxy 2 Section High-Low Back Section	8
Accessories	9
Warranty	10





Face hole on head section



Fold-down armrests

Two-layer padding



Powder coated, extremely stable scissor steel frame

Easy access foot switch

Paper Roll Holder

High-performance, multidirectional casters

Personal Authorization System



Easy access foot switch





Galaxy 5 Section High-Low



Product features:

- Fold-down armrests in head section
- Adjustable head section from +40 to -70 with gas spring
- Easy access foot switch, accessible from any position around the table 360°, for simple height adjustment
- Face hole on head section with removable face hole cover
- Adjustable leg section from 0 to 90 with gas spring
- Two-layer padding with non-flammable, bio-compatible and scratch resistant upholstery in many colours
- Powder coated, extremely stable scissor steel frame
- 2 side brackets for fixation of stabilization belts and elastic bands
- Traction support sleeves
- Personal Authorization System with 2 safety magnetic keys
- Hallotronic® System electric high/low height adjustment from 52 to 102 cm
- Non-slip rubber feet with 1 cm height adjustment allowing leveling of the table
- 2 types of frame:
 - 2 small wheels and 2 legs
 - 4 directional castors with central brake and 4 legs

Table Articulations/Dimensions

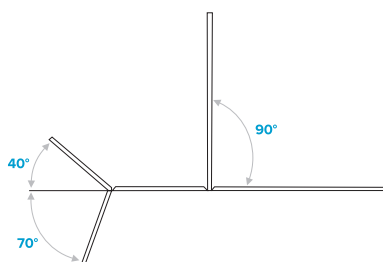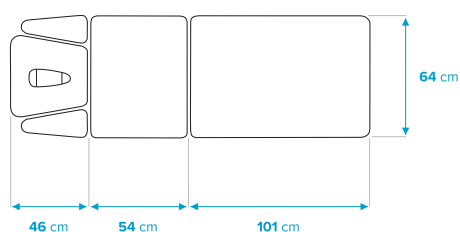

Ordering Information

Part Number Description

Galaxy 5 Section High-Low Treatment Table

3153-X	Galaxy 5 Section Electrical, All-round Switch, Castors
3152-X	Galaxy 5 Section Electrical, All-round Switch
3151-X	Galaxy 5 Section Hydraulic, Castors
3150-X	Galaxy 5 Section Hydraulic

Technical Specifications

Electrical

Mains Power:	100 -240 Volt 50/60 Hz
Standby Consumption:	max. 1 A
Weight:	Standard 92 kg/4 Castors 104 kg
Shipping Weight:	Table Weight + 50 kg
Shipping Dimensions*:	220 cm x 85 cm x 77 cm
Lifting Capacity:	200 kg
Height Adjustment:	52 - 102 cm
Electrical safety class:	Class II, Type B

Hydraulic Pump

Weight:	Standard 93 kg/4 Castors 102 kg
Shipping Weight:	Table Weight + 50 kg
Shipping Dimensions*:	220 cm x 85 cm x 77 cm
Lifting Capacity:	200 kg
Height Adjustment:	51 - 96 cm

Upholstery Color Options

		
Ecru-Crème (x = 00)	Beige (x = 01)	Blue (x = 02)
		
Navy Blue (x = 03)	Terra Cotta (x = 04)	Chocolate (x = 05)
		
Gray (x = 06)	Graphite Gray (x = 07)	Black (x = 08)

NOTICE: Color charts are produced within the limitations of lithographic color printing and are intended as a general guide only.

For ordering purposes, please use model number, hyphen, upholstery color number as indicated in the chart above.

Galaxy 4 Section High-Low

Ordering Information	
Part Number	Description
Galaxy 4 Section High-Low Treatment Table	
3143-X	Galaxy 4 Section Head-Rest, Electrical, All-round Switch, Castors
3142-X	Galaxy 4 Section Head-Rest, Electrical, All-round Switch
3141-X	Galaxy 4 Section Head-Rest, Hydraulic, Castors
3140-X	Galaxy 4 Section Head-Rest, Hydraulic

Technical Specifications	
Electrical	
Mains Power:	100 -240 Volt 50/60 Hz
Standby Consumption:	max. 1 A
Weight:	Standard 86 kg/4 Castors 98 kg
Shipping Weight:	Table Weight + 50 kg
Shipping Dimensions*:	220 cm x 85 cm x 77 cm
Lifting Capacity:	200 kg
Height Adjustment:	52 - 102 cm
Electrical safety class:	Class II, Type B
Hydraulic Pump	
Weight:	Standard 87 kg/4 Castors 96 kg
Shipping Weight:	Table Weight + 50 kg
Shipping Dimensions*:	220 cm x 85 cm x 77 cm
Lifting Capacity:	200 kg
Height Adjustment:	51 - 96 cm

Upholstery Color Options		
		
Ecru-Crème (x = 00)	Beige (x = 01)	Blue (x = 02)
		
Navy Blue (x = 03)	Terra Cotta (x = 04)	Chocolate (x = 05)
		
Gray (x = 06)	Graphite Gray (x = 07)	Black (x = 08)


NOTICE: Color charts are produced within the limitations of lithographic color printing and are intended as a general guide only. For ordering purposes, please use model number, hyphen, upholstery color number as indicated in the chart above.



Product features:

- Fold-down armrests in head section
- Adjustable head section from +40 to -70 with gas spring
- Easy access foot switch, accessible from any position around the table 360°, for simple height adjustment
- Face hole on head section with removable face hole cover
- Two-layer padding with non-flammable, bio-compatible and scratchproof upholstery in many colours
- Powder coated, stable scissor steel frame
- 2 side brackets for fixation of stabilization belts and elastic bands
- Traction support sleeves
- Personal Authorization System with 2 safety magnetic keys
- Hallotronic® System electric high/low height adjustment from 52 to 102 cm with capacity up to 200 kg
- Non-slip rubber feet with 1 cm height adjustment allowing levelling of the table
- 2 types of frame:
 - 2 small wheels and 2 legs
 - 4 directional castors with central brake and 4 legs

Table Articulations/Dimensions

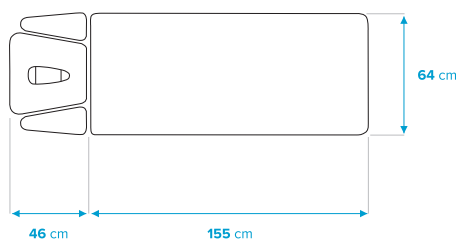



Galaxy 3 Section High-Low



Product features:

- Adjustable head section from +40 to -70 with gas spring
- Easy access foot switch, accessible from any position around the table 360°, for simple height adjustment
- Face hole on head section with removable face hole cover
- Adjustable leg section from 0 to 90 with gas spring
- Two-layer padding with non-flammable, bio-compatible and scratch resistant upholstery in many colours
- Powder coated, extremely stable scissor steel frame
- 2 side brackets for fixation of stabilization belts and elastic bands
- Traction support sleeves
- Personal Authorization System with 2 safety magnetic keys
- Hallotronic® System electric high/low height adjustment from 52 to 102 cm
- Non-slip rubber feet with 1 cm height adjustment allowing leveling of the table
- 2 types of frame:
 - 2 small wheels and 2 legs
 - 4 directional castors with central brake and 4 legs

Table Articulations/Dimensions

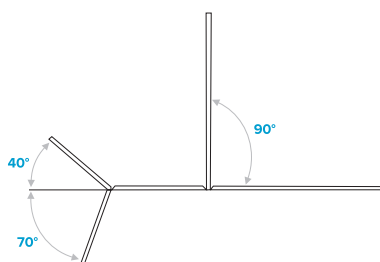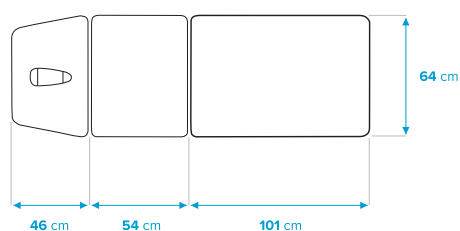

Ordering Information

Part Number Description

Galaxy 3 Section High-Low Treatment Table

3133-X	Galaxy 3 Section Electrical, All-round Switch, Castors
3132-X	Galaxy 3 Section Electrical, All-round Switch
3131-X	Galaxy 3 Section Hydraulic, Castors
3130-X	Galaxy 3 Section Hydraulic

Technical Specifications

Electrical

Mains Power:	100 - 240 Volt 50/60 Hz
Standby Consumption:	max. 1 A
Weight:	Standard 87 kg/4 Castors 99 kg
Shipping Weight:	Table Weight + 50 kg
Shipping Dimensions*:	220 cm x 85 cm x 77 cm
Lifting Capacity:	200 kg
Height Adjustment:	52 - 102 cm
Electrical safety class:	Class II, Type B

Hydraulic Pump

Weight:	Standard 88 kg/4 Castors 97 kg
Shipping Weight:	Table Weight + 50 kg
Shipping Dimensions*:	220 cm x 85 cm x 77 cm
Lifting Capacity:	200 kg
Height Adjustment:	51 - 96 cm

Upholstery Color Options

		
Ecru-Crème (x = 00)	Beige (x = 01)	Blue (x = 02)
		
Navy Blue (x = 03)	Terra Cotta (x = 04)	Chocolate (x = 05)
		
Gray (x = 06)	Graphite Gray (x = 07)	Black (x = 08)

NOTICE: Color charts are produced within the limitations of lithographic color printing and are intended as a general guide only.

For ordering purposes, please use model number, hyphen, upholstery color number as indicated in the chart above.

Ordering Information

Part Number Description

Galaxy 2 Section High-Low Treatment Table

3123-X	Galaxy 2 Section Head-Rest, Electrical, All-round Switch, Castors
3122-X	Galaxy 2 Section Head-Rest, Electrical, All-round Switch
3121-X	Galaxy 2 Section Head-Rest, Hydraulic, Casters
3120-X	Galaxy 2 Section Head-Rest, Hydraulic

Technical Specifications

Electrical

Mains Power:	100 -240 Volt 50/60 Hz
Standby Consumption:	max. 1 A
Weight:	Standard 81 kg/4 Castors 93 kg
Shipping Weight:	Table Weight + 50 kg
Shipping Dimensions*:	220 cm x 85 cm x 77 cm
Lifting Capacity:	200 kg
Height Adjustment:	52 - 102 cm
Electrical safety class:	Class II, Type B

Hydraulic Pump

Weight:	Standard 82 kg/4 Castors 91 kg
Shipping Weight:	Table Weight + 50 kg
Shipping Dimensions*:	220 cm x 85 cm x 77 cm
Lifting Capacity:	200 kg
Height Adjustment:	51 - 96 cm

Upholstery Color Options

		
Ecru-Crème (x = 00)	Beige (x = 01)	Blue (x = 02)
		
Navy Blue (x = 03)	Terra Cotta (x = 04)	Chocolate (x = 05)
		
Gray (x = 06)	Graphite Gray (x = 07)	Black (x = 08)

NOTICE: Color charts are produced within the limitations of lithographic color printing and are intended as a general guide only. For ordering purposes, please use model number, hyphen, upholstery color number as indicated in the chart above.


Galaxy 2 Section High-Low Head



Product features:

- Adjustable head section from +40 to -70 with gas spring
- Face hole on head section with removable face hole cover
- Easy access foot switch, accessible from any position around the table 360°, for simple height adjustment
- Two-layer padding with non-flammable, bio-compatible and scratch resistant upholstery in many colours
- Powder coated, stable scissor steel frame
- 2 side brackets for fixation of stabilization belts and elastic bands
- Traction support sleeves
- Personal Authorization System with 2 safety magnetic keys
- Hallotronic® System electric high/low height adjustment from 52 to 102 cm with capacity up to 200 kg
- Non-slip rubber feet with 1 cm height adjustment allowing leveling of the table
- 2 types of frame:
 - 2 small wheels and 2 legs
 - 4 directional castors with central brake and 4 legs

Table Articulations/Dimensions

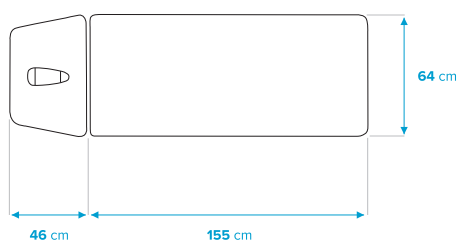



Galaxy 2 Section High-Low Back



Product features:

- 2 section massage and treatment table
- Adjustable back section from +?? to -??with gas spring
- Face hole on head section with removable face hole cover
- Easy access foot switch, accessible from any position around the table 360°, for simple height adjustment
- Two-layer padding with non-flammable, bio-compatible and scratch resistant upholstery in many colours
- Powder coated, stable scissor steel frame
- 2 side brackets for fixation of stabilization belts and elastic bands
- Traction support sleeves
- Personal Authorization System with 2 safety magnetic keys
- Hallotronic® System electric high/low height adjustment from 52 to 102 cm with capacity up to 200 kg
- Non-slip rubber feet with 1 cm height adjustment allowing levelling of the table
- 2 types of frame:
 - 2 small wheels and 2 legs
 - 4 directional castors with central brake and 4 legs

Table Articulations/Dimensions

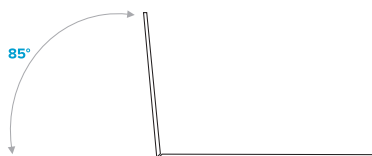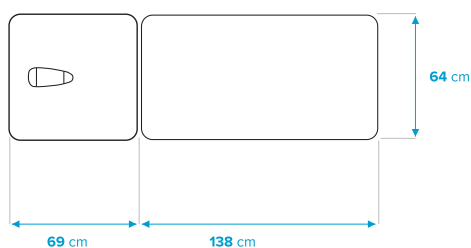

Ordering Information

Part Number Description

Galaxy 2 Section High-Low Back Treatment Table

3127-X	Galaxy 2 Section Back-Rest, Electrical, All-round Switch, Castors
3126-X	Galaxy 2 Section Back-Rest, Electrical, All-round Switch
3125-X	Galaxy 2 Section Back-Rest, Hydraulic, Castors
3124-X	Galaxy 2 Section Back-Rest, Hydraulic

Technical Specifications

Electrical

Mains Power:	100 -240 Volt 50/60 Hz
Standby Consumption:	max. 1 A
Weight:	Standard 93 kg/ 4 Castors 105 kg
Shipping Weight:	Table Weight + 50 kg
Shipping Dimensions*:	220 cm x 85 cm x 77 cm
Lifting Capacity:	200 kg
Height Adjustment:	52 - 102 cm
Electrical safety class:	Class II, Type B

Hydraulic Pump

Weight:	Standard 94 kg/ 4 Castors 103 kg
Shipping Weight:	Table Weight + 50 kg
Shipping Dimensions*:	220 cm x 85 cm x 77 cm
Lifting Capacity:	200 fkg
Height Adjustment:	51 - 96 cm

Upholstery Color Options

		
Ecru-Crème (x = 00)	Beige (x = 01)	Blue (x = 02)
		
Navy Blue (x = 03)	Terra Cotta (x = 04)	Chocolate (x = 05)
		
Gray (x = 06)	Graphite Gray (x = 07)	Black (x = 08)

NOTICE: Color charts are produced within the limitations of lithographic color printing and are intended as a general guide only.



For ordering purposes, please use model number, hyphen, upholstery color number as indicated in the chart above.

Ordering Information	
Part Number	Description
Accessories	
3110GALXYPPRHLDR	(A) Paper Roll Holder for the Galaxy Tables
3711-xx	(B) Traction Supports (2 pack)
3712	(C) Stabilization Belt
3713	(D) Kaltenborn Wedge
3714-xx	(E1) Wedge
3715-xx	(E2) Wedge
3716-xx	(F1) Bolster
3717-xx	(F2) Bolster
3718-xx	(F3) Bolster

Upholstery Color Options		
		
Ecrú-Crème (x = 00)	Beige (x = 01)	Blue (x = 02)
		
Navy Blue (x = 03)	Terra Cotta (x = 04)	Chocolate (x = 05)
		
Gray (x = 06)	Graphite Gray (x = 07)	Black (x = 08)

NOTICE: Color charts are produced within the limitations of lithographic color printing and are intended as a general guide only. For ordering purposes, please use model number, hyphen, upholstery color number as indicated in the chart above.

Accessories

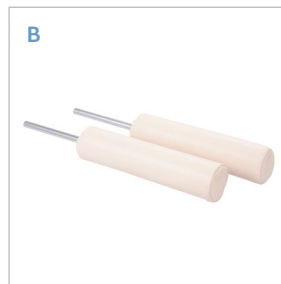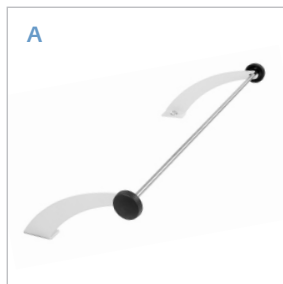

Product features:

- A** Paper Roll Holder - Holds paper rolls up to 60 cm wide
- B** Traction Supports - Mechanical traction supports
- C** Stabilization Belt - Available in black, this stabilization belt attaches to side brackets (available on most models)
- D** Kaltenborn Wedge - For use by professionals trained in manual mobilization techniques.
- E1 & E2** Wedge - Wedges and Bolsters are used to reinforce the correct positioning of the patient during treatment.
- F1, F2 & F3** Bolster

Warranty

DJO, LLC ("Company"), warrants that the Galaxy Tables ("Products") are free of defects in material and workmanship. This warranty shall remain in effect for two years (24 months) for the cushions and three years (36 months) for the base frame and motor from the date of original consumer purchase. If this Product fails to function during the warranty period due to a defect in material or workmanship, at the Company's option, the Company or selling dealer will repair or replace this Product without charge within a period of thirty days from the date on which the Product is returned to the Company or the dealer.

All repairs to the Product must be performed by a service center certified by the Company. Any modifications or repairs performed by unauthorized centers or groups will void this warranty.

This Warranty Does Not Cover:

- Replacement parts or labor furnished by anyone other than the Company, the selling dealer, or a service technician certified by the Company.
- Defects or damage caused by labor furnished by someone other than Company, the selling dealer, or a certified Company service technician.
- Any malfunction or failure in the Product caused by product misuse, including, but not limited to, the failure to provide reasonable and required maintenance or any use that is inconsistent with the Product User Manual.

DJO, LLC is not responsible for injury or damage resulting from modifications or service performed by non-authorized DJO, LLC service personnel.

COMPANY SHALL NOT BE LIABLE IN ANY EVENT FOR INCIDENTAL OR CONSEQUENTIAL DAMAGES.

Some locations do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation or exclusion may not apply to you.

To obtain service from Company or the selling dealer under this warranty:

1. A written claim must be made within the warranty period to the Company or the selling dealer.
2. The Product must be returned to the Company or the selling dealer by the owner.

This warranty gives you specific, legal rights and you may also have other rights which vary from location to location.

The Company does not authorize any person or representative to create for it any other obligation or liability in connection with the sale of the Product.

Any representation or agreement not contained in the warranty shall be void and of no effect.

THE FOREGOING WARRANTY IS IN LIEU OF ALL OTHER WARRANTIES, EXPRESSED OR IMPLIED, INCLUDING ANY WARRANTY OR MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE.

Galaxy Tables are manufactured by Meden-Inmed Sp. z o.o. exclusively for Chattanooga. A DJO Global Company.



DJO GLOBAL

AUSTRALIA:

T: +1300 66 77 30
F: +1300 66 77 40
E: customerservice.au@DJOglobal.com

BENELUX:

T: Belgium 0800 18 246
T: Netherlands 0800 0229442
T: Luxemburg 8002 27 42
E: benelux.orders@DJOglobal.com

CANADA:

T: +1 1866 866 5031
F: +1 1866 866 5032
E: canada.orders@DJOglobal.com

CHINA:

T: (8621) 6031 9989
F: (8621) 6031 9709
E: information_china@DJOglobal.com

DENMARK, FINLAND, NORWAY & SWEDEN:

T: Sweden 040 39 40 00
T: Norway 8006 1052
T: Finland 0800 114 582
T: Denmark +46 40 39 40 00
E: info.nordic@DJOglobal.com

FRANCE:

T: +33 (0)5 59 52 86 90
F: +33 (0)5 59 52 86 91
E: sce.cial@DJOglobal.com

GERMANY:

T: +49 761 4566 01
F: +49 761 456655 01
E: infoservice@DJOglobal.com

ITALY:

T: +39 02 484 63386
F: +39 02 484 09217
E: vendite@DJOglobal.com

INDIA:

T: +91 44 66915127
E: customercare.india@DJOglobal.com

SOUTH AFRICA:

T: +27 (0) 87 3102480
F: +27 (0) 86 6098891
E: info.southafrica@DJOglobal.com

SPAIN:

T: +34 934 803 202
F: +34 934 733 667
E: ventas@DJOglobal.com

SWITZERLAND:

T: +41 (0) 21 695 2360
F: +41 (0) 21 695 2361
E: info@compex.ch

UK & IRELAND:

T: +44 (0)1483 459 659
F: +44 (0)1483 459 470
E: ukorders@DJOglobal.com

UNITED STATES:

T: +1 800 336 6569
F: +1 800 936 6569
E: customercare@DJOglobal.com

DJO GLOBAL, EXPORT CENTRES

ASIA-PACIFIC:

DJO Asia-Pacific Limited
Unit 1905, 19/F, Tower II
Grand Central Plaza
138 Shatin Rural Committee Road
Shatin
HONG KONG
T: +852 3105 2237
F: +852 3105 1444
E: info.asia@DJOglobal.com

EUROPE, MIDDLE EAST & AFRICA:

DJO Benelux
Welvaartstraat 8
2200 Herentals
BELGIUM
T: +32 (0) 14248350
F: +32 (0) 14248358
E: info.emea@DJOglobal.com

LATIN AMERICA:

DJO Global, Inc
1430 Decision Street
Vista
CA 92081-8553
U.S.A.
T: 1 800 336 6569
F: 1 800 936 6569
E: info.latam@DJOglobal.com



Together in Motion™