



Genesis Walker

The lightest full shell walker boot on the market



Genesis Walker Tall & Mid-Calf



Genesis Walker

Not all walker boots are created equal.

Inspired by **athletic shoe design technology**, the Genesis Walker provides a **biomechanical, natural walk**. The result is a walker boot that mimics a **normal shoe and healthy gait**, allowing you to recover with minimal disruption to your daily life.



The lightweight, low profile **Genesis Walker** revolutionizes patient **comfort, maneuverability and control**.



The Genesis Advantage

Comfort

Maneuverability

Comfort with Every Step

The **Genesis Walker** is the **lightest full shell pneumatic walker on the market***. It has been computer engineered to provide strength where required while eliminating extra weight.



Dynamic Air Stride Technology
Similar to popular athletic shoes, the Genesis Walker utilizes a **dynamic air chamber** in the heel of the tread which actively **absorbs the shock of heel strike** and allows for a more comfortable walking experience.

Trim-able Shell

The Genesis Walker may be adjusted for patients with larger calf muscles. A cut line allows the boot to be cut to accommodate a large calf.



Happy feet, Healthy body

Narrow Heel Width

The wider the heel width, the more force that can be transferred to the knee joint. The Breg heel is designed to **mimic the width of a normal shoe**. This makes it **easier for the patient to walk** in the boot because it is similar to the shoe on the opposite foot.



15°

Proper heel height and width reduces pelvic tilt and secondary knee pain

Low Heel

The Genesis Walker has the **lowest heel height on the market** making it comparable to a standard shoe**. Uneven heel height can cause pelvic tilt resulting in secondary pain in the hip and back. By having a comparable heel to a standard shoe, the Genesis Walker reduces pelvic tilt which **helps reduce secondary hip and back pain**.

* Data on file.

** Standard shoe is defined as a standard men's dress shoe with an approximate heel height of 1 inch.

Control


Walk Normally While Recovering

The Genesis Walker's **streamlined design** provides a low profile look and feel. This reduces the feeling of bulk and makes it noticeably easier to maneuver.

Engineered Tread Design

Rocker Bottom

Breg studied the components of a normal gait. The tread **mimics the action of the ankle** during the gait cycle. Because in a walker boot your ankle is immobilized, this rocker bottom tread mimics the actions of a healthy ankle **promoting a more natural walk**.

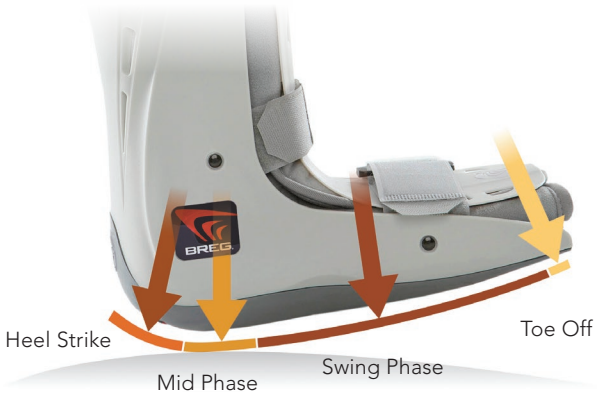


Scalloped No-Slip Tread
Athletic shoe tread design provides traction on a variety of surfaces permitting better cross-directional maneuverability.

Performance that Delivers

The Genesis Walker delivers control, support and protection to the foot and ankle from time of injury until full recovery.

Dual Air Chambers

Two identical air chambers provide compression from the foot to ankle while maintaining ankle position within the boot. This allows for increased control and comfort while healing occurs.

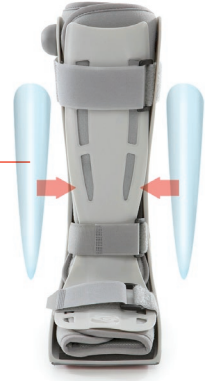

Built-in Pneumatic System

Built-in Air Pump

The built-in air pump allows patient to make incremental increases or decreases in pressure around the foot and ankle. The built-in pump system eliminates the need to carry a hand pump.



Durable, Hard Shell Design

The hard shell design combines the control of a cast with the convenience of a walker. It also allows adjustments for swelling and inspection of the foot and ankle. The hard shell design provides enhanced protection to an injury while healing. The boot can also be removed for foot and ankle inspection or dressing changes.



Genesis Mid-Calf Walker

Utilizing the same great features as the Genesis Walker, the Genesis Mid-Calf Walker is designed for maximum **comfort**, **manueverability** and **control**, but in a more compact size.



Great things come in **small packages**

Using the same cutting-edge design as the tall Genesis Walker, the mid-calf boot's short stature allows it to provide **optimal support** with **minimal impact** to your life.



Mid-Calf Design

At 25% shorter than the standard Genesis Walker, the Genesis Mid-Calf is low profile and ultra lightweight.

The lightest full shell walker boot just got lighter



Available with **custom logos**

Genesis Walker Tall

Genesis Walker Mid-Calf

Ordering Information

Customer Care

At Breg, we are committed to total customer care. Each of our friendly Customer Care Representatives receives extensive training on all Breg products and practices, ensuring that you always get the best service possible. From our diverse line of medical devices to our exceptional staff, we strive to exceed your expectations.

To order, call: 800-321-0607 | 800-897-BREG(2734)

Fax: 800-329-2734 | 800-959-2734

Available Monday - Friday, 6:00 am to 5:00 pm, Pacific Time

International

International customers can order by calling: +1-760-795-5440

Part Numbers

Part Number	Size	Shoe Size		Product
BL515001	Extra Small	M: 1-2.5	W: 1-3.5	Genesis Walker 3-Strap
BL515003	Small	M: 2.5-5.5	W: 3.5-6	Genesis Walker 3-Strap
BL515005	Medium	M: 5.5-9	W: 6-9.5	Genesis Walker 3-Strap
BL515007	Large	M: 9-12.5	W: 9.5-13	Genesis Walker 3-Strap
BL515009	Extra Large	M: 12.5-17	W: 13+	Genesis Walker 3-Strap
BL514001	Extra Small	M: 1-2.5	W: 1-3.5	Genesis Walker 4-Strap
BL514003	Small	M: 2.5-5.5	W: 3.5-6	Genesis Walker 4-Strap
BL514005	Medium	M: 5.5-9	W: 6-9.5	Genesis Walker 4-Strap
BL514007	Large	M: 9-12.5	W: 9.5-13	Genesis Walker 4-Strap
BL514009	Extra Large	M: 12.5-17	W: 13+	Genesis Walker 4-Strap
BL525001	Extra Small	M: 1-2.5	W: 1-3.5	Genesis Mid-Calf Walker
BL525003	Small	M: 2.5-5.5	W: 3.5-6	Genesis Mid-Calf Walker
BL525005	Medium	M: 5.5-9	W: 6-9.5	Genesis Mid-Calf Walker
BL525007	Large	M: 9-12.5	W: 9.5-13	Genesis Mid-Calf Walker
BL525009	Extra Large	M: 12.5-17	W: 13+	Genesis Mid-Calf Walker



To learn more about Breg,
visit www.breg.com

2885 Loker Ave. E
Carlsbad, CA 92010
(800) 897-BREG(2734)